



アデノ随伴ウイルス(AAV)  
Adeno-associated virus

U220713B

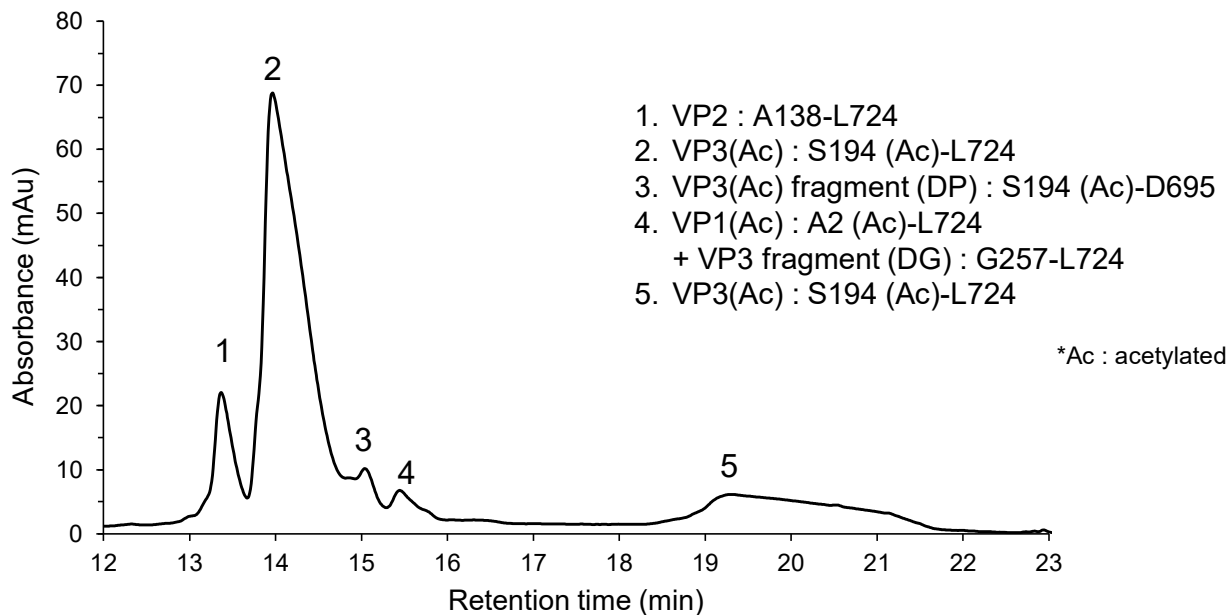

**Sample preparation**

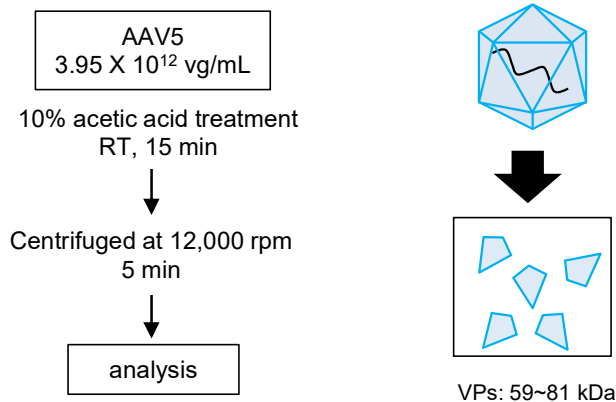

Courtesy of Prof. S. Uchiyama, Osaka University  
This research was supported by AMED under Grant Number JP18ae0201001.

|             |   |
|-------------|---|
| Column      | : Accura Triart Bio C18 (1.9 μm, 30 nm)<br>150 X 2.1 mmI.D.   |
| Eluent      | : A) water/difluoroacetic acid (100/0.1)<br>B) acetonitrile/difluoroacetic acid (100/0.1)<br>20-32%B (0-1 min), 32-36%B (1-16 min), 36-80%B (16-20 min) |
| Flow rate   | : 0.2 mL/min  |
| Temperature | : 80°C  |
| Detection   | : UV at 280 nm<br>ESI-MS (positive ion mode)  |
| Injection   | : 50 μL   |
| Sample      | : Denatured AAV5  |